



ALL IN ONE MEDIA PREPARATION



SIMPLIFY

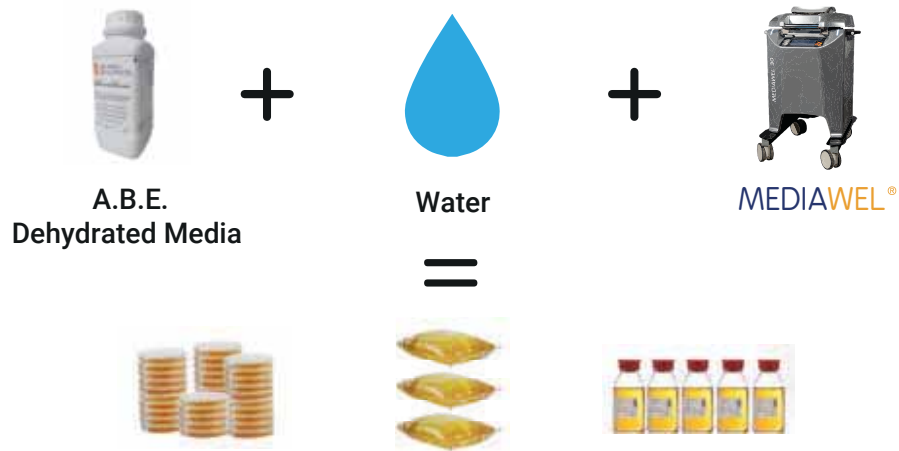


HOMOGENEITY

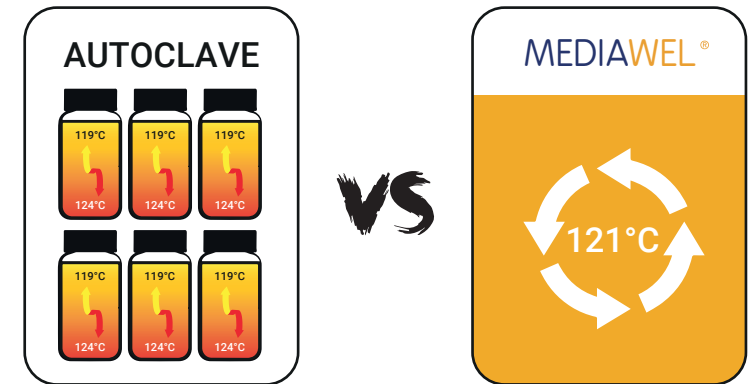


ALL IN ONE MEDIA PREPARATION

TIME SAVING

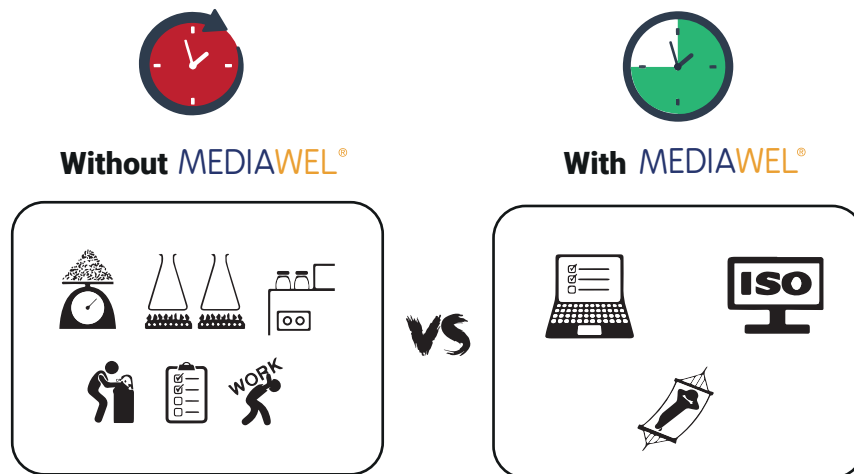

FAST COOLING SYSTEM

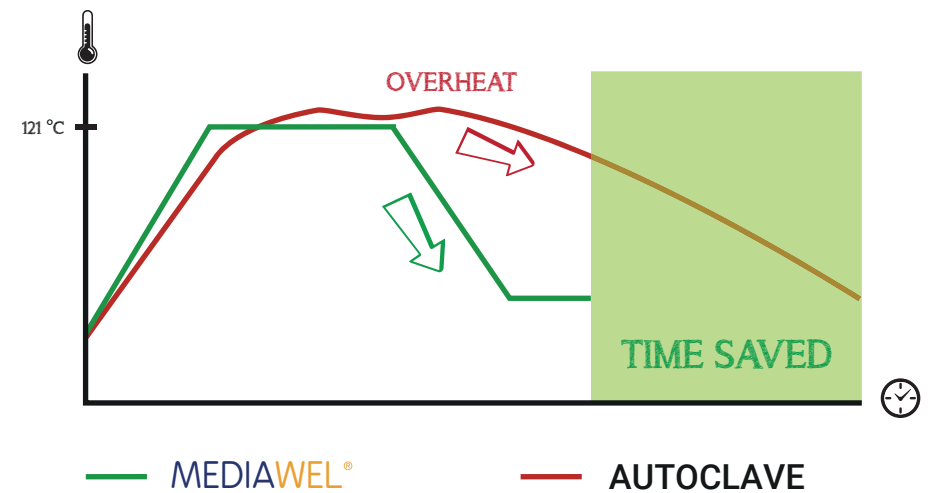

MEDIA PREPARATOR

All in one media preparation !

MEDIAWEL 10[®]

Double magnetic stirring blade :

Large surface for an optimum temperature homogeneity

Standardize your media

With 50 programs available



Transparent safety lid

Locked above 80°C

Audit trail

CFR21 PART 11

Bactystopper[®] filtration

Inside the vessel



e-TRACE[®] technology



7" LCD touchscreen

Gloves friendly



Multicolored visual alarm

Data export by USB or in real time through the Labpage[®] online service



Up to 10 000 cycles saved

Trolley with medical wheels

Lockable on the front, easily slid under the bench



Easy cleaning

No valve, removable bucket

Direct and smooth evacuation with dedicated space for a bowl



PERFORMANCES



From 1 to 12 liters



45 - 75 min



25-123 °C / ± 0.5 °C



TRACEABILITY WITH YOUR MEDIAWEL®

Each cycle is exported (PDF/CSV) to your USB key :



REAL TIME DATA MONITORING Labpage® Real time data monitoring through WIFI/Ethernet, connected to your PC / Smartphone / Tablet :

Dashboard for the live monitoring : program settings, temperature settings and current temperature, alarms and notifications ...



Cycles are automatically recorded and stored. Recall them one by one or sort them out by batch number, date, user ...



**SD card included - up to 10 000 cycles data storage*

Delivered with the MEDIAWEL 30 : User Manual, Control certificate, Large magnetic stirrer (MEDW2008), Sampling tube (MEDW2009), Draining plug (MEDW2010), Sampling nozzle seals x2 (MEDW4011), Sampling nozzles x2 (MEDW2003 : Ø 12mm + MEDW2004: Ø 6.4mm), Inlet and outlet water tubing (4 m) + 4 tubing clamps (MEDW4012), Decompression filters x2 (MEDW3001), Bactystopper - Calibrated sterile carded cotton x50 (MEDW3003), Sliding washer stirrer axis x2 (MEDW4003), Sampling nozzle seals x2 (MEDW4011), USB key, 3 Months free access to the Labpage "Real time data monitoring LICENCE"

TECHNICAL DATA SHEET

Capacity	From 2.5 to 30 liters
Sterilization temperature	From 25 to 125°C
Sterilization time	From 1 to 99 min
Distribution temperature	From 25 to 80°C
Temperature accuracy	+/- 0.5°C after stabilization
Remaining volume	<130ml
Addition cap	✓
Traceability	eTRACE® Module integrated, USB data export (PDF/CSV)
Traceability WIFI/Ethernet (optional)	Labpage® dashboard. Real time data monitoring, data storage on SD card (up to 10 000 cycles)
Temperature probe	Media : 3 probes (1 standard, 1 in option and 1 for calibration)
Safety pressure release	Safety valve
Dispensing port	Adapters for tubes of 6.4 mm / 12 mm
Continuous homogenization	✓
Touchscreen	Tactile, LCD 7"
Alarm	Multicolored visual alarm
Cleaning	Trapdoor, easy access to the vessel
High mobility	Stable and easy to move : 4 large medical wheels
Robustness	Antishock body (ABS PC PLASTIC)
Vessel	Stainless steel AISI 316L
Programs	50
Weight	154kg
Dimensions (L x l x H)	735 x 585 x 1053 mm
Warranty	1 year
Power supply	400V +/- 10% - 50/60 Hz - 17A (by phase) - 9000 W or 208V +/- 10% - 50/60 Hz - 9A (by phase) - 2070 W
Country of origin	France - European Union


MEDIAWEL 30®

Quick water connection kit (MEDW2002) :




Quick and easy water connection from one room to another !

MEDIAWEL Dust cover (MEDW2021) :

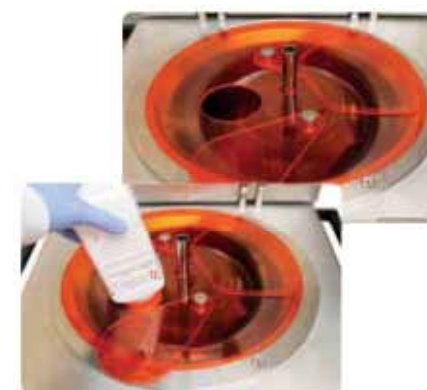

Avoid volatile powder !


FOOD SAMPLE PREPARATION

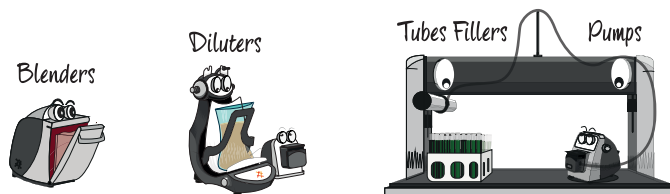

AUTOMATED MEDIA PREPARATION




SAMPLE COLLECTION

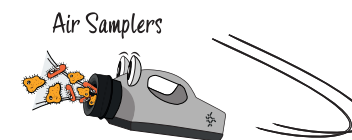

BOOK YOUR DEMO !

<https://www.alliance-bio-expertise.com/en/store-locator/>

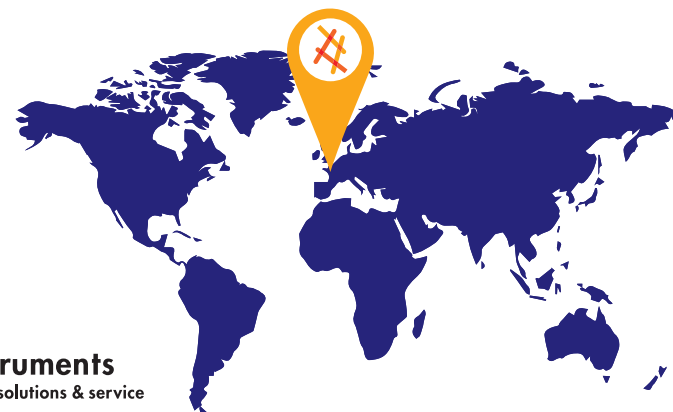

Instruments
smart solutions & service

IGZ Instruments AG
Furtbachstrasse 17
8107 Buchs ZH

Tel. +41 44 456 33 33
igz.ch igz@igz.ch



FOLLOW THE BLOG

[alliance-bio-expertise.com](https://www.alliance-bio-expertise.com)



ZI de Courbouton, Le Tremplin - 35480 GUIPRY - FRANCE

E-mail : contact@abioexpertise.com

Phone : +33(0)2.40.51.79.53

[alliance-bio-expertise.com](https://www.alliance-bio-expertise.com)